



CLASS II BIOSAFETY CABINET BEV1J7

CLASS II BIOSAFETY CABINET BEV1J7

BIOLOGICAL SAFETY CABINET

Biological safety cabinets (BSC) is clean equipment in biological experiments, to prevent damage of the experimental operation of poisonous and harmful airborne particles, in order to protect the safety of the operator and the environment in the experiment. BSC provides reliable security. Standard with a HEPA filter to protect products and environment, plenum leak prevention technology, innovation stable inflow speed system, with LCD screen, microcomputer control technology, especially for biological contamination of a variety of high-risk experiments.

Meet the domestic YY0569-2005 BSC requirements and NSF / ANSI 49-2008, the European Union EN12469: 2000 BSC standard, to provide users with a satisfactory safety assurance.

IQ, OQ, PQ and other related certification services.

CE certification. Two-year warranty period.



Clean Air Technology

Equipped with German EBM high airflow and low energy consumption fan, based on airflow mechanics principle, arc-shaped large-diameter inlet design to reduce turbulence generated, to ensure low noise operation of the system.

Reasonable duct design, optimization of the plenum tank, with the design of micro-shaped mesh layer outflow wind deflector, non-convex, no dead corners to ensure uniformity and stability of the airflow.

Equipped with long life exhaust HEPA and downflow ULPA, >99.995% efficient at 0.12 μm particles.

BSC have two types: A2 Type: 70% inner circulation, 30% exhaust; B2 Type: 100% exhaust. Work zone width: 1000 mm and 1300 mm.

Low energy consumption motor

SPECIFICATIONS

Model	BEV1J7
Clean rating	International standard ISO14644.1 CLASS 5 (US209E, CLASS 100)
Inflow velocity	0.56 m/s
Down-flow velocity	0.35 m/s
Down-flow velocity precision	± 0.015 m/s
Noise level (dB)	≤ 65
Vibration half-peak value	$\leq 3 \mu\text{m}$ (X, Y, Z direction)
High efficient filter	with ULPA super efficient filter, efficiency 99.999%;
Light intensity (LX)	≥ 650
Pressure difference range (Pa)	0-500
Working environment	Environment temp. 10-30°C, relative humidity under 70%, no obvious vibration and dust
External dimensions (HxWxD, mm)	2210x1120x850
Working zone dimension (HxWxD, mm)	700x1000x500
Power supply voltage	Single phase AC220V/50HZ
Power (W)	300
Net weight (KG)	175
Gross weight (KG)	195
Package size (HxWxD, mm)	1720x1290x1110
Alt Name	Biological Safety Cabinet

ACCESSORIES FOR PURCHASE

No	Name
1	Exhaust pipe and fan(no installation accessories)

FEATURES

Clean Air Technology

Equipped with German EBM high airflow and low energy consumption fan, based on airflow mechanics principle, arc-shaped large-diameter inlet design to reduce turbulence generated, to ensure low noise operation of the system.


Reasonable duct design, optimization of the plenum tank, with the design of micro-shaped mesh layer outflow wind deflector, non-convex, no dead corners to ensure uniformity and stability of the airflow.

Equipped with long life exhaust HEPA and downflow ULPA, >99.995% efficient at 0.12 μm particles.


BSC have two types: A2 Type: 70% inner circulation, 30% exhaust; B2 Type: 100% exhaust. Work zone width: 1000 mm and 1300 mm.



Low energy consumption motor

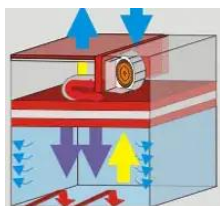


A2 (70%inner circulating, 30%discharge)

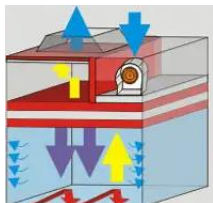

B2 (100%discharge)

Stable Fan Speed System

High air flow fan, the system will automatically track and adjust the air volume and pressure, to balance the role of compensation, with the design of micro-shaped mesh layer outflow wind deflector, so airflow has more stability and uniformity.

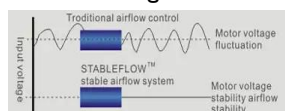

stable airflow system




Micro-mesh lamina flow board

Plenum Leak-proof Technology

Around the work zone at the top and side wall plenum design, so that the discharge of polluting gases completely in a plenum environment, equipped with indoor backplane to a molding design, is an effective way to avoid the possibility of leakage.

Equipped with the accurate pressure sensor to monitor the pressure of the working area, where there is abnormal situation such as breakage of the filter zone or negative pressure seal failure, the system will automatically stop the fan running, and through sound and light alarm prompts the operator, thus avoiding due to leakage caused damage.

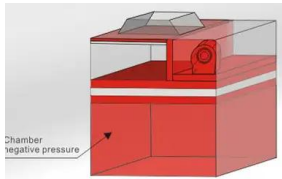


Plenum leak-proof technology

Microcomputer Control

The integrated LCD screen. Display inflow/downflow speed, pressure, filter life, failure alarm and other parameter status.

The user-friendly design, with sterilization, lighting interlock function, front window overhigh alarm and front window closed then the fan shut down.

Fan speed is adjustable according to the need to adjust the appropriate amount of fan.

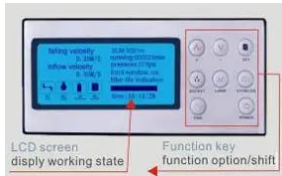

Excellent Imported Sensor: Air Velocity Sensor And Pressure Difference Sensor

Imported air velocity sensors collect data of high precision, the performance of corrosion-resistant, anti-interference ability, real-time monitoring and control cabinet status.

Imported pressure difference sensor is responsive, effectively monitoring the pressure of the work zone and plenum.



Pressure difference sensor



Air velocity sensor





Biolab Scientific Ltd.

Trillium Executive Center, East Tower, 675 Cochrane Dr, Markham, Ontario L3R 0B8, Canada
Email: info@biolabscientific.com | Website: www.biolabscientific.com