



-86°C ULTRA LOW TEMPERATURE UPRIGHT FREEZER BJX-202

-86°C ULTRA LOW TEMPERATURE FREEZER BJX-202

-86°C ULT FREEZER



Temperature Control System

Microprocessor control, Large touch screen display inner temperature clearly, and with easy viewing.

The inner temperature can be adjustable at the range of 40°C~-86°C, with an increment of 0.1°C.

Lock and password for protection to prevent arbitrary adjustment of operating parameters.

Optional 7" color LCD touch screen display system.

Safety Protection Measures

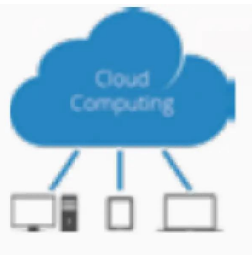
Malfunction alarms: High-temperature alarm, low-temperature alarm, sensor failure alarm, power failure alarm, low backup-battery voltage alarm, and over-temperature alarm. Alarm temperatures can be set as required.

SPECIFICATIONS

Model	BJX-202
Capacity	280L
Temperature range	-40~-86°C
Cooling performance(°C)	-80
Rated Voltage	220/50
Exterior Material	Galvanized steel powder coating
Interior Material	Galvanized steel powder coating
Insulation Material	PUF+VIP
Climate Class	N
Cooling Type	Direct cooling
Defrost Mode	Manual
Refrigerant	HC mixed
Cabinet Type	Upright
Power (W)	520
Energy consumption (KW.H/24H)	9.5
Shelves (stainless steel)	3
Caster optional	Yes
Test Hole	Yes
Temperature Recorder	Optional
Door locking device	Yes
Handle	Yes
Pressure balance hole	Yes
Racks & Boxes	Optional
Display	7" Touch screen with USB port
Low/High temperature alarm	Yes
Condenser high temperature alarm	Yes
Power failure alarm	Yes
Sensor failure alarm	Yes
High room temperature alarm	Yes
Alarm mode	Sound and light alarm, remote alarm terminal
Thickness of Cabinet Foamed Layer	90mm

Inner Door Number	2
Power Supply(V/Hz)	220/50
Capacity for 2 inch boxes	180
Net/Gross Weight(approx)	190kg/213kg
Interior Dimensions(WDH)	460x470x1310mm
Exterior Dimensions(WDH)	715x725x1980mm
Package Dimensions(WDH)	735x750x2190mm
Thickness of Door	90mm
Alt Name	-86°C ULT Freezer

ACCESSORIES FOR PURCHASE

No	Image	Name
1		CryoRacks
2		CryoBoxes
3		Co2/LN backup device
4		Chart temperature recorder
5		Wireless Monitoring system

FEATURES

Temperature Control System

Microprocessor control, Large touch screen display inner temperature clearly, and with easy viewing.
The inner temperature can be adjustable at the range of 40°C~86°C, with an increment of 0.1°C.
Lock and password for protection to prevent arbitrary adjustment of operating parameters.
Optional 7" color LCD touch screen display system.

Safety Protection Measures

Malfunction alarms: High-temperature alarm, low-temperature alarm, sensor failure alarm, power failure alarm, low backup-battery voltage alarm, and over-temperature alarm. Alarm temperatures can be set as required.
Three alarm methods: Buzzer alarm, digital flashing alarm, and remote alarm.
Protection functions: Power-on delay and interval protection between stops.

Refrigeration System

Optimized cascade refrigeration technology with a SECOP compressor for effective cooling.
Specially designed, high-efficiency heat exchanger for improved heat transfer and energy savings.
VIP insulation panels ensure excellent thermal insulation performance.
CFC-free, overall polyurethane foam insulation.
High-efficiency hydrocarbon refrigerant.

Ergonomic Design

Safety door lock, preventing unauthorized access.
Wide-voltage design from 192 V to 242 V.
Double inner doors reduce cold loss after opening the main door.
Suitable for 2-inch and 3-inch international-standard cryo boxes.
Pressure-balanced port design for easy door opening.
Heavy-duty casters with adjustable locking screws for easy movement and secure positioning.
Ultra-simple condenser-filter design for convenient maintenance and cleaning.

APPLICATIONS

-86°C ULT Freezer is specially designed for long term storage of various biological products, such as germs, virus, erythrocytes, leukocytes, cutis. It can be installed in institutions including blood banks, hospitals, epidemic prevention services, research institutes and laboratories for electronic and chemical plants, biological engineering institutes and marine fishery companies. From individual laboratory users to ever-increasing sample storage needs, you can choose from ultra-low temperature storage solutions in multiple sizes.



Biolab Scientific Ltd.

Trillium Executive Center, East Tower, 675 Cochrane Dr, Markham, Ontario L3R 0B8, Canada
Email: info@biolabscientific.com | Website: www.biolabscientific.com