



## NANO SPECTROPHOTOMETER BSNA-102

# NANO SPECTROPHOTOMETER BSNA-102

## MICRO-SPECTROPHOTOMETER



### User-Friendly Software, Easy to Use

Graphical software operation, more intuitive interface, the results can be directly exported, easy to save, view and output data.

### Micro-Volumes Measuring

Only 0.5  $\mu$ L-2  $\mu$ L sample is needed for each test. After the measurement, the samples can be recovered and the precious samples can be studied with confidence.

### Long Life Light Source, Do Not Need to Warm Up

Xenon flash, life span is  $10^9$  (up to 10 years). No preheating, direct use, ready to test at any time.

### High Concentration Detection

The maximum concentration of the detectable sample is 15000 ng/ $\mu$ L, and the sample basically does not need to be diluted.

## SPECIFICATIONS

Model	BSNA-102
Wavelength range	200~800 nm
Minimum sample size	0.5~2.0 $\mu$ L
Path length	0.2 mm 1.0 mm
Light source	Xenon flash lamp
Detector type	2048-linear CCD array
Wavelength accuracy	1 nm
Spectral resolution	$\leq 3$ nm
Absorbance precision	0.003 Abs
Absorbance accuracy	1 % (7.332 Abs at 260 nm)
Absorbance range	0.04~90 A
Nucleic acid detection range	2~4500 ng/ $\mu$ L (dsDNA)
Measurement time	< 5 s
Dimension (WxDxH) mm	200x250x166
Weight	2.6 kg
Sample pedestal material	Aluminum alloy and quartz fiber
Operating voltage	DC 24 V 2 A
Operating power	20 W
Standby power	5 W
Software compatibility	Win 7, Win XP, Win 8
Cuvette Mode (OD600 Measurement)	
- Light source	/
- Wavelength range	/
- Absorbance range	/
Fluorometer Mode	
- Sensitivity	/
- Linear dynamic range	/
- Repeatability	/
Alt Name	Micro-Spectrophotometer

## FEATURES

### User-Friendly Software, Easy to Use

Graphical software operation, more intuitive interface, the results can be directly exported, easy to save, view and output data.

### Micro-Volumes Measuring

Only 0.5  $\mu\text{L}$ -2  $\mu\text{L}$  sample is needed for each test. After the measurement, the samples can be recovered and the precious samples can be studied with confidence.

### Long Life Light Source, Do Not Need to Warm Up

Xenon flash, life span is  $10^9$  (up to 10 years). No preheating, direct use, ready to test at any time.

### High Concentration Detection

The maximum concentration of the detectable sample is 15000 ng/ $\mu\text{L}$ , and the sample basically does not need to be diluted.

### Convenient and Easy to Use


Directly point the sample on to the sample plate without dilution or cuvette. The sample concentration can be measured as 50 times of the conventional uv-visible photometer, and the result can be directly output as the sample concentration.

### Fast Detection

No dilution or cuvette needed in the detection process; 5 s can complete the test and display the result.

### Single Machine Operation, Convenient and Efficient

## APPLICATIONS

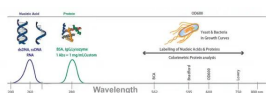

### Different Fluorescence Channels:

#### UV Channel:

Excitation wavelength: 365 $\pm$ 20 nm

Common reagent: Hoechst 33258, 4-MU, EnzChek Caspase

Nucleic acid quantification, plant GUS reporter gene detection, apoptosis detection

#### Blue Channel:

Excitation wavelength: 460 $\pm$ 20 nm

Common reagent: PicoGreen, oligreen, RiboGreen, GFP, Protein, Fluorescein

dsDNA, ssDNA, GFP, gene detection, fluorescein detection, protein quantification.

#### Green Channel:

Excitation wavelength: 525 $\pm$ 20 nm

Common reagent: Rhodamine, Cy3, RFP Vybrant Cytotoxicity

Rhodamine detection, Cy-3 fluorescence labeling detection, RFP gene detection, cytotoxicity detection.

#### Red Channel:

Excitation wavelength: 625 $\pm$ 20 nm

Common reagent: Cy5, Quant-iT RNA

Cy-5 fluorescence labeling detection, RNA quantification

## MORE INFO

It is a micro-spectrophotometer with full wavelength (200-800 nm).

The detectable sample is 50 times that of conventional UV / Vis spectrophotometer, and maximum detection concentration up to 4500 ng/ $\mu\text{L}$  (dsDNA).

High resolution CCD array detectors, 5 s can complete the detection, display the results.

Long life pulse xenon lamp light source.

The need to connect PC computers to run detection, data saving, printing, sorting are very convenient.



**Biolab Scientific Ltd.**

Trillium Executive Center, East Tower, 675 Cochrane Dr, Markham, Ontario L3R 0B8, Canada  
Email: [info@biolabscientific.com](mailto:info@biolabscientific.com) | Website: [www.biolabscientific.com](http://www.biolabscientific.com)